B-Digital Baby Scale

(Manual - Handleiding - Mode d'emploi - Gebrauchsanleitung- Instrucciones de uso - Instruções de utilização - Istruzioni per l'uso - Návod k použití)

B400800

Made In China – Bo Jungle - Chamo BV – Brambroek 23B - 9500 GERAARDSBERGEN – BELGIUM

www.bojungle.eu – office@bojungle.eu – Tel: +32 54 240 331 – Fax: +32 9 270 01 80

! WARNING

IMPORTANT! PLEASE READ CAREFULLY & KEEP FOR FUTURE REFERENCE

Thank you for choosing our product. Please take a little time to read the important safety notes detailed below, as this will ensure many years of happy, safe use.

This electronic scale is a device for checking the baby's weight, it is especially useful for weighing newborns and premature babies, where the weight is an exceptionally important development factor. We used our best efforts to make sure the product is on top quality, comfortable and safe during weighing.

Device instruction

1. ON OFF button
 - Press the button to switch on or off the device
 - If you use a towel to lay the baby on, put the towel on the platform and press the button to set the display to "0.00"
2. UNIT button
 - Press the button to adjust the unit of measurement : gram(g),kilogram(kg), pound(lb)or stone(st:lb)
3. Display (LCD)
 - The size: 85x 28mm

Product specification

- 1) Device dimensions: 630X300X54mm
- 2) The Capacity: 50g ~ 20KG
- 3) Press ON, press OFF/auto OFF
- 4) Accuracy: d=0.005KG
- 5) 4 high precision strain gauge sensor
- 6) Tare function
- 7) Low battery and over load indicate
- 8) Auto lock/hold function (easier weighing of moving babies-the result is automatically locked weight)
- 9) Large LCD: 85x28mm
- 10) Battery: 3 x 1.5v AAA battery (not included)

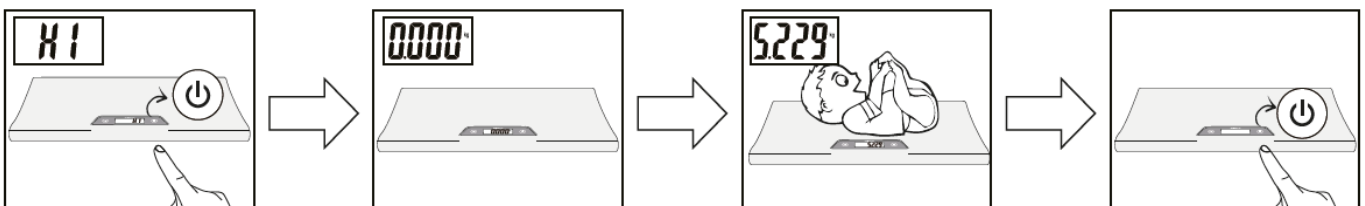
Battery installation

- 1) Before using the scale, check the battery compartment and remove any plastic insulation strips.
- 2) Insert the batteries (3 x 1.5V AAA) into the battery compartment.
- 3) If you are replacing the batteries,make sure the positive and negative contacts are properly aligned.

Device operation

Weigh




- 1) Please use the scale on the hard & flat surface. Do not weight on carpet or any soft surface as it may affect accuracy.
- 2) Press "⏻" button, the display shows "Hi" ;
- 3) After display reads "0.00kg" then it could be able to measure weight ;
- 4) Put baby or other stuff on the scale platform, and LCD display will auto lock the weight after 2 seconds;
- 5) Press "⏻" key to turn off the scale , or scale will turn off itself after 30 seconds.

Reset and tare function



It could weight continuously:

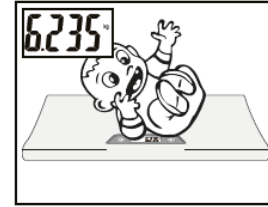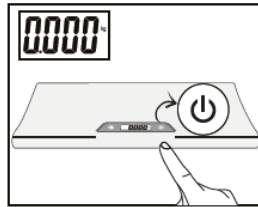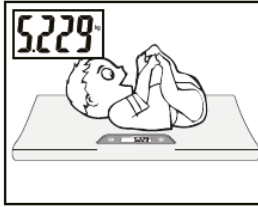- 1) First, put the first one baby on the scale platform to get weight;
- 2) Second, press "⏻" button, display will turn back to zero, then put the next one on the scale. (this function is repeatable).

Unit conversion

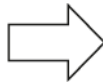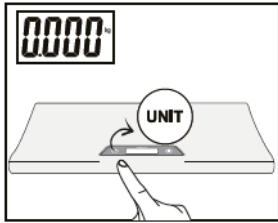Press "UNIT" button to convert "kg" "g" "lb" and "st:lb".

Overload Indication

If the weight is over loaded "max capacity" LCD display will show "ERR".

Low Power Indication

When LCD display shows "LO", Please replace the new battery with correct polarity.

Safety precautions

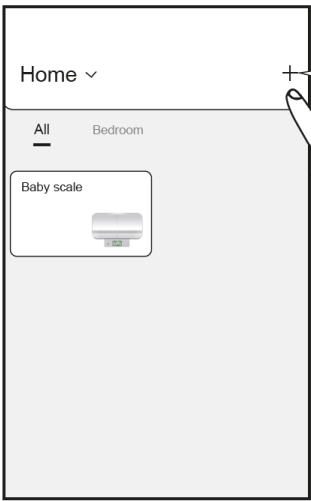
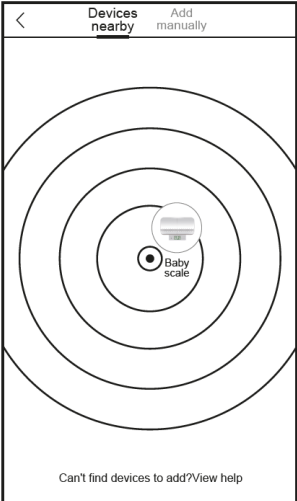
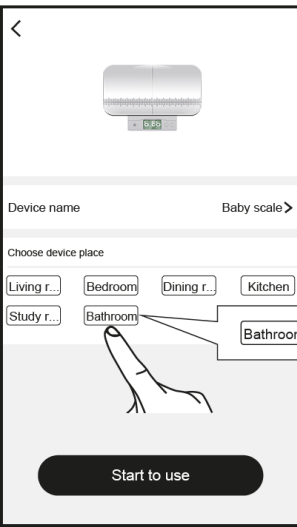
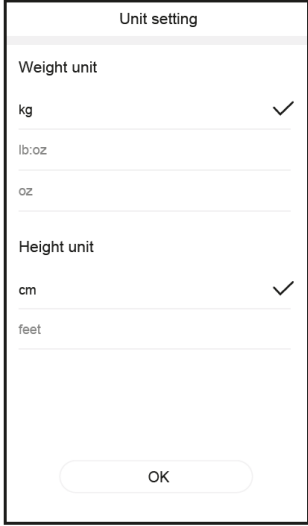
- When in use, the scale should be placed on flat, dry and hard surface.
- Use to this device for its intended purpose at home.
- Do not overload the scale's capacity, the max.is 20kgs.
- Avoid shaking or hitting the scale.
- Do not submerge in water.
- When in storage, do not place anything on the scales
- The optimal temperature for operating and storing the scales is from 5C to 30C.

Download App

APP name: AiLink

Remark: Our APP is available for IOS and android device with IOS version 8.0 or Android version 5.0 and Bluetooth 4.0 (Or higher version)	Download APP "AiLink" on Apple store or Google play online store

Connect with AiLink APP

<p>1. Open the app and the bluetooth of the phone. After entering the app, click the “+” button in the upper right corner as shown in figure 1.</p>		<p>2. The APP automatically scans nearby devices when entering the interface to add devices, as shown in figure 2.</p>	
<p>3. Click the scanned device to add automatically, and it will be displayed after successful addition as shown in figure 3.</p>		<p>4. Select the weighing position on the app and click “start”.</p>	
<p>5. Set the units on the APP.</p>			

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